

SAFETY LANCETS CATALOG



Key Features / Benefits: Lightweight

Small and lightweight design reducing production and shipping costs which are passed on to our customers.

Safe

Engineered to prevent accidental Finger-sticks with two redundant safety features:

- 1) an automatically retracting needle that cannot be reused.
- 2) a release button that remains depressed after use to ensure one-time operation.

Comfortable

Split-second incision by beveled needle minimizes pain and helps increase patient comfort.

Easy to Use


Intuitively designed and ready to use, the lancet is easily operated and activated by a gentle push of the release button.

Directions for Use:

1. Clean puncture site.
 2. Remove protective cap from the Safety Lancet.
 3. Place the lancet against the puncture site.
 4. Push the release button.
 5. Dispose lancet after use.
- Do not use if the seal has been damaged or broken.

Code	Color	Diameter	Depth
4101	●○ mazarine+white	33G	1.5 mm
4102	●○ Pink+ wathet blue	30G	1.5 mm
4103	●○ wathet blue+white	28G	1.5 mm
4104	●○ wathet blue+white	27G	1.5 mm
4105	●○ Pink+Purple	26G	1.8mm
4106	●○ Pink+white	25G	1.8mm
4107	●○ Green+white	23G	1.8mm
4108	●○ Orange+white	21G	1.8mm


Package : 20Boxes/Carton

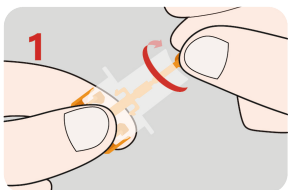

1. Remove protective cap from the Safety Lancet.



3. Place the lancet against the puncture site.

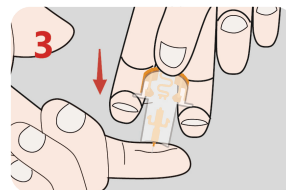

4. Push the release button.



5. Dispose lancet after use. Do not use if the seal has been damaged or broken.